


DISPENSING PACIFIER

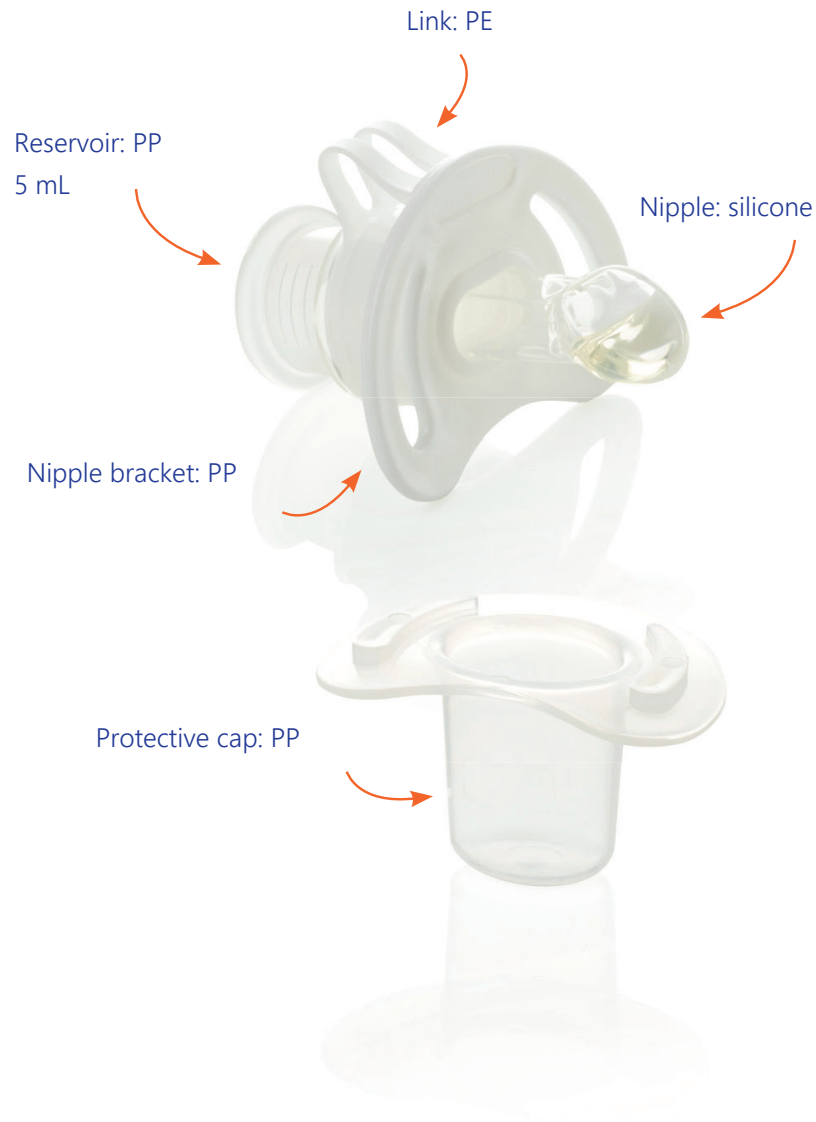

With the Dispensing Pacifier, administering liquid is no longer a fight.

The dispensing pacifier is a comfortable way to give young patients any prescriptions they need, with no added stress for parents. Perfectly adapted for pediatric use, this device makes it a lot easier to administer liquid supplements and treatments. The pacifier consists of five pieces: a nipple with a reservoir and a protective cap, a bracket and a link.

USE

1. Unscrew the reservoir and fill it with the liquid according to the prescribed dose
2. Screw the reservoir on the nipple bracket
3. Remove the protective cover and put the pacifier into the baby's mouth